

The no-drilling solution: Icon® Caries Infiltrant

Icon is a breakthrough, minimally invasive way your dentist can stop the progression of an early cavity before it needs to be treated using a needle or drill. It fills and seals the tiny holes and stabilizes the tooth, offering many unique advantages:

- The treatment does not require drilling the tooth or shots of anesthesia
- Healthy tooth structure is preserved
- The treatment is completed in just one visit to your dentist
- The treated surface of the tooth looks like healthy enamel – no more white spots!
- By treating cavities early, the life expectancy of the tooth can be prolonged

How a cavity becomes a cavity...

What causes you to get a cavity? Bacteria that live in your mouth thrive on sugars and starches, producing acids as a result. These acids eat away at the enamel on your teeth, creating tiny holes or pores.

This eventually weakens the structure of the tooth, often appearing as a chalky white spot on the surface of the tooth. As the process progresses, a cavity is formed that your dentist traditionally treats by drilling and filling the tooth.

Are you a good candidate for Icon® treatment?

Icon is the breakthrough treatment that can...

- Stop an early cavity from progressing before it needs to be “drilled and filled”
- Cosmetically remove white spots caused by early cavities or orthodontic treatment

Keep in mind, Icon is only effective for treating cavities in early stages. Cavities that have progressed beyond the early stage will need to have a conventional filling or other treatment.

Ask your dentist if you are a good candidate for treatment with Icon!

drilling-no-thanks.info
[#IconByDMG](https://www.instagram.com/IconByDMG)



Patients prefer Icon®

A new way to treat early cavities and remove white spots without drilling

- ✓ No anesthesia
- ✓ No drilling
- ✓ One visit



Icon Proximal


Quick and easy treatment for early-stage cavities!

Icon penetrates into the tiny holes in tooth enamel created by acids in your mouth. It then fills and reinforces the enamel to block the cavity from progressing further. The cavity is prevented from progressing any further! Here's how the treatment works...

Early stage cavities between the teeth are detected using an x-ray or other diagnostic tool. To treat these with Icon, your dentist will first use a tool to separate the teeth slightly. The affected tooth will be gently and painlessly prepared with a gel. The dentist will then use an ultra-thin applicator to fill the early cavity with Icon.



How it works

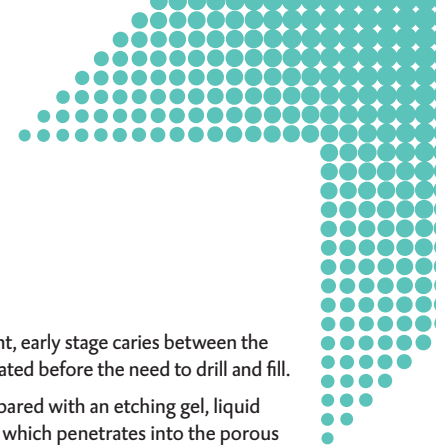


Icon Proximal



1. Before treatment, early stage caries between the teeth can be treated before the need to drill and fill.
2. After being prepared with an etching gel, liquid resin is applied, which penetrates into the porous tooth tissue where it cures under light.
3. The porous tissue is filled and blocks the incipient caries from progressing

Icon Smooth Surface

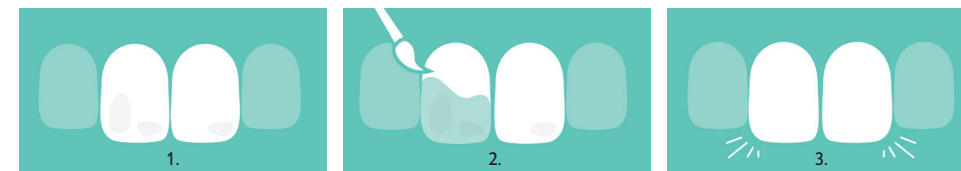

1. Before treatment, caries, fluorosis or enamel developmental defects have caused unesthetic discoloration.
2. After being prepared with an etching gel, liquid resin is applied, which penetrates into the porous tooth tissue where it cures under light.
3. The porous tissue is filled. The infiltrated area adapts optically to the healthy tooth thanks to a similar refractive index of infiltrant and enamel.

Icon Smooth Surface

Removal of white spots in one visit, with no drilling!

Early-stage cavities on the smooth surface of a tooth are usually detected visually. They often appear as chalky white spots that can alter the appearance of your teeth. These white spots also commonly appear on teeth after orthodontic treatment.

Icon provides a fast and easy cosmetic treatment for these type of white spots. The affected tooth is prepared with a gel. The dentist will then fill the affected area with Icon. White spots treated with Icon take on the appearance of the surrounding healthy enamel, restoring the natural look of the tooth.

No drilling or anesthesia is needed to cosmetically remove white spots with Icon. And, unlike other treatments, you only need to make one visit to the dentist!



White spots before Icon treatment.



After Icon treatment, the beautiful, natural look of the treated teeth are restored!